




# WASHING DISINFECTION DEVICE ACCESSORIES



[www.sumerinc.com](http://www.sumerinc.com)

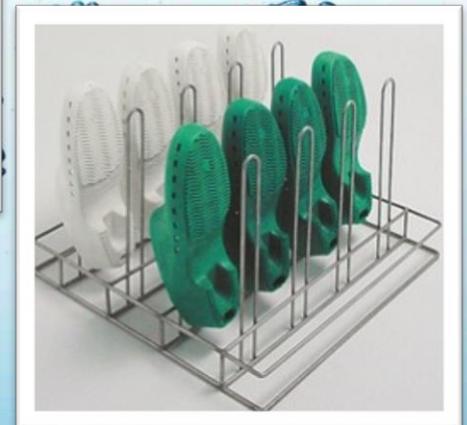



## About Us

Sümer A.Ş. was established in 1981 in Ankara to provide services in the medical device sector. It has aimed advancement since the day of its establishment by also taking growth and compliance with the contemporary technologies and protecting the environmental conditions. Our manufacturing takes place in an area of 25.000 m<sup>2</sup> in Ankara Organized Industrial Zone.

Our Firm is following the innovations in its sector and in abroad through its research and development unit and its application staff with a strong infrastructure of engineers, and is continuing to produce devices it had developed in computer environment based on such innovations with high technology and to contribute their development so as to be most beneficial for the Turkish medicine.

With this purpose, our Firm is strictly following the “Quality Management” principles and rules from design of the products to the after-sale servicing.

Our Firm has been currently certificated for compliance with ISO 9001 quality management system, ISO 13485 medical device quality management system certificate and ISO 14001 environment management system certificate and with product certificates under MDD 93/42/EEC Medical Devices Directive and PED 2014/68/EUPressurized Equipment certificate. Furthermore, our steam sterilizers, and disinfection and washing devices have been certificated by the German accredited body.

Our Firm is successfully implementing several projects supported by Tübitak, KOSGEB, and Ankara Development Agency. SÜMER A.Ş. has been awarded with the following:  
Ankara Industry Chamber Award for 2012, 2011 R&D Award,  
2012 Success Award.

Our Firm possesses the following certificates;  
CE Certificates under the following directives:  
ISO 9001,  
ISO 13485,  
ISO 14001 Quality Management System,  
MDD 93/42/EEC Medical Devices,  
PED 2014/68/EUPressurized Equipment



## Our Vision

To make the Sümer brand a global brand to make our Firm remembered first in the sector.

## Our Mission

Our main task is to create designs with competition power in the global sense by taking the priorities of the sector into consideration and being respectful to the environment and people and giving the first priority to the wishes and expectations of customers, and also to produce innovative technological medical products by meeting all the national and international legal requirements.

## Our Basic Values


- We are bound up with the Medical Ethical Rules,
- We are people-oriented,
- We respect environment,
- We are creative,
- We are customer-oriented,
- We are innovative,
- We are pro-active,
- We believe in the team spirit.

## WASHING DISINFECTION DEVICE ACCESSORIES

**ADELA** series washing disinfection devices offer a comprehensive range of devices and accessories for automatic washing and disinfection. Our automated surgical cleaning experience in Central Sterile Feeding Departments, endoscopy centers and special applications enables safe, efficient and cost-effective cleaning and disinfection procedures based on proven modern technologies.



## Washer Disinfector Accessories

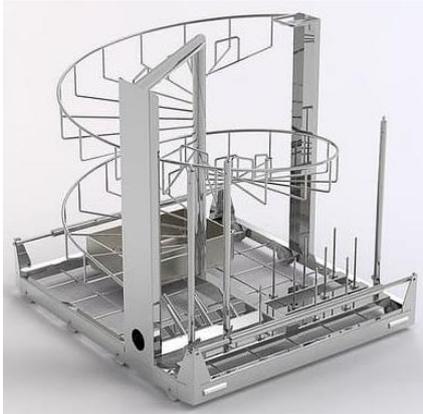

Anaesthesia Rack

Ref. No: 8698918589343


AN spiral wash carts offer a high-capacity solution with its 12 hose connections for hoses of an inner connection  $\varnothing$  of 8-30 mm and a maximal length 1500 mm. An additional 15 injectors (L=102 mm), 15 injectors (L= 30 mm), 6 injectors for breathing bags and 2 injectors for bellows are mounted. There is also an insert on the central pipe for mounting a tray for masks



Instrumental Rack

Ref. No: 8698918589148

This wash cart is designed to maximize productivity in processing standard instrument trays with 3,4,5 shelves for loading (four removable shelves and spray arms).

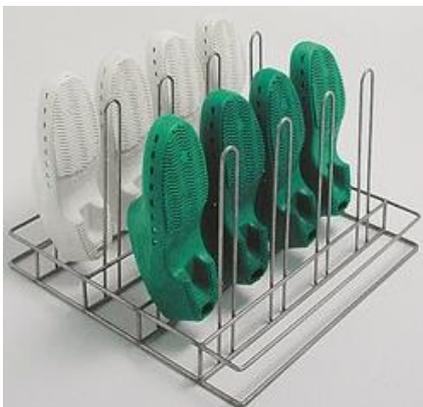

Shoes Rack

Ref. No: 8698918589356

This rack holds up to 15 pairs of shoes for washing and disinfection at the same time. It is placed in an instrument wash cart.



Container Rack

Ref. No: 8698918589543

The container cart is designed to assure efficient cleaning and disinfection. Furthermore, the inclined positioning of the goods ensures maximum run-off and short drying times.

## Washer Disinfector Accessories



Baby Bottle Rack

Ref. No: 869891855325

This stand for baby bottles is designed to optimize capacity while delivering great cleaning and disinfection results. The stand is equipped with handles for easy loading and unloading of the cart.



MIS Rack

Ref. No: 86989122147

For limited throughput of tubular instruments, we recommend an I-module. This could be well suited in a TSSU or medical clinic.



Tray Rack

Ref. No: 8698918589509

The tray rack holds up to 9 trays.



Emesis Basin Tray

Ref. No: 8698918589313

The tray rack holds up to 9 emesis basin trays.




## Dental Washer Accessories



Loading Basket

Ref. No: 869891853532

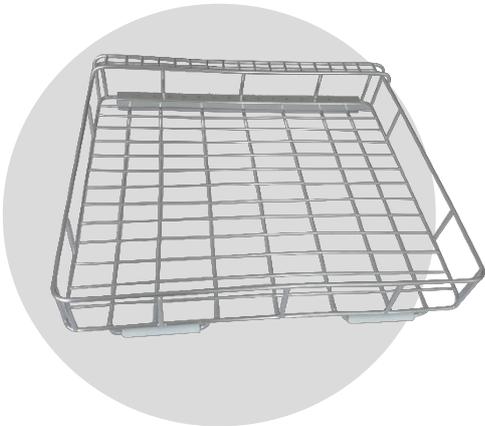


Basic Washing Basket

Ref. No: 86989148202



Vertical Basket For Mirrors And Instruments

Ref. No: 8698945563

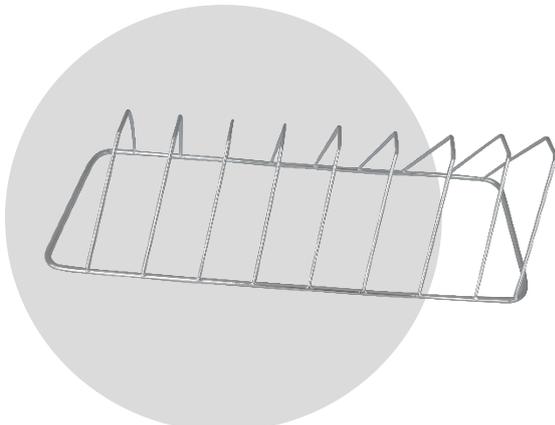

8 Positions Tray Basket

Ref. No: 8698958903



#### Solution Partners

SIEMENS

Schneider  
Electric

seko

HANNING

OSV

SKF

DIGIEMSA

WÜRTH

SICK



[www.sumerinc.com](http://www.sumerinc.com)

 **sumer**

Address : Başkent OSB Başkent Bulvarı No:81  
Malıköy-Sincan/Ankara/Türkiye  
Phone : +90 312 394 61 73 +90 312 503 08 21 (pbx)  
Fax : +90 312 394 61 95  
E-mail : [sumer@sumeras.com](mailto:sumer@sumeras.com)

 **SZUTEST**  
NOTIFIED BODY 2195

 **kiwa**  
partner for progress

 CE

2195

 ISO  
13485

 ISO  
9001

 ISO  
14001